



SMART SKIN TECHNOLOGY.

genio
SMART XS

Operating instructions

genio
SMART SKIN TECHNOLOGY.


Thank you very much for deciding for genio - SMART SKIN TECHNOLOGY.

The unique combination of ultrasound and red/blue light therapy, combined with carefully developed sophisticated serums, provides professional facial skincare therapy at home!

Ten minutes every second day is all it takes to get healthy, youthful-looking skin!

Please read these instructions carefully for optimum results and maximum treatment success.

Overview

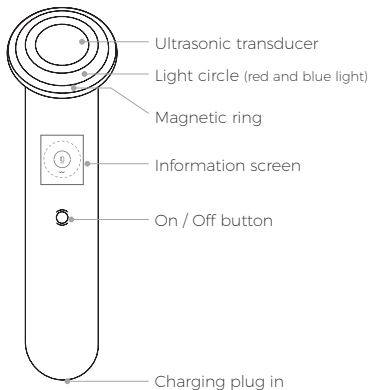

Introduction

1. Short press the On / Off button to enter the nutrition mode. The program starts with a smooth vibration and the ultrasound waves and the red / blue lights work simultaneously. To optimize the light therapy, the colour changes every minute.
2. Press the button again to shut down.
3. When the battery power is low (only 10 minutes or less) the outer blue circle starts blinking. Please recharge.
4. When charging: the outer blue circle flashes. After being fully charged: the outer blue circle stops blinking.
5. After 5 minutes the vibrating mode will be more intensive to show that you have already reached 50% of the treatment time.
6. The countdown (from 10 to 1) on the information screen shows how many minutes of the treatment time are still left.
7. After the 10 minutes the device will automatically shut down.


Display



The program

Ultrasonic waves combined with red / blue light therapy based on different wave lengths and 3 MHz power, provides the optimized physical interaction for best results.

Accessories content SMART XS

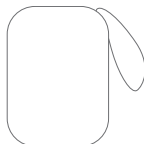

Travel bag



Device



Magnetic ring



Charging cable



Manual


The genio therapy



1. Cleanse the face as you do in your daily routine.



2. Spray 2-3 pumps of the facial water in the area you want to treat.

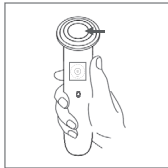

3. Apply the serum on the transducer.




4. Switch on the device by pressing the On / Off button.

The genio therapy



4. Move the device in slow motions over your face and neck. The program lasts for 10 min (optimal treatment time for face and neck). Slow movements will support the absorption of the agents.

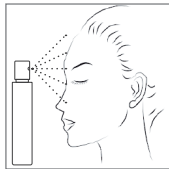

5. When the serums gets dry (not sliding smoothly anymore), spray a few pumps of the facial water again in the area which you are treating. This will optimize the effect and support the Ultrasound effect and treatment success.



6. After the operation, the remaining serum can be massaged in with your fingers. You can use 1-2 pumps of the facial water to support this process.

TIP: You can use the genio therapy to see sustainable results, every second day. The program will lift up your skin, diminish big pores, work against inflammation, and remove fine lines and wrinkles. You will recognize a visible change in your skin after 2 weeks.

Please do not extend the treatment time of 10 minutes per day/per treatment in one area (f.e face).

Product specification

Product name	genio SMART XS
Product model	BM21
Product weight	118G
Product size	181x52x44mm
Battery capacity	700mAh
Product power	3W
Accessories	Magnetic ring, cable, manual, travel bag

Frequently asked questions and answers

Please check our website: genio-innovations.com for further information.

Environmental instructions

1. This product meets the requirements of "Management Methods for Restriction of the use of Hazardous Substances in Electrical and Electronic Products".
2. Within the environmental protection period, it will not have leakage and precipitation of harmful substances that will affect consumers' health during normal use. Please use it with confidence.
3. The environmental protection use period of our products is 10 years, and the "environmental protection use period" is only valid when the product is used under the normal conditions described in this manual.
4. In order to protect the environment, please separate the batteries and other parts of this product from domestic garbage and send them to a qualified recycling station for dismantling. Dismantle and recycle the product according to relevant national regulations.

Environmental instructions

The name and content of harmful substances in products						
Part name/hazardous substance	lead	HG	cadmium	Hexavalent chromium	PBBs	PBDE
Plastic components	○	○	○	○	○	○
Circuit board (electronic) components	×	○	○	×	○	○
Connecting wire	○	○	○	○	×	×
Silicone Rubber	○	○	○	○	○	○
Screws and washers	○	○	○	×	○	○
Motor	○	○	○	○	○	○
Battery	○	○	○	○	○	○

This table is prepared in accordance with the provisions of SJ/T 11364.

○: indicates that the content of the hazardous substance in all homogeneous materials of the part is below the limit requirement stipulated in GB/T26572.

×: indicates that the content of the hazardous substance in at least one homogeneous material of the part exceeds the limit specified in GB/T26572.

Warning

- To ensure safe and correct use of this appliance, do not operate the unit without thoroughly reading this manual first.
- While using implantable medical electronic devices such as heart pacemakers and artificial heart and lung life-supporting electronic devices is prohibited to use this device.
- Patients with fever, infectious diseases and malignant tumors are prohibited to use.
- Do not stay on the same part for more than 1 minute to avoid skin problems.
- This product has a built-in lithium battery, please do not put it in fire, and avoid using in high temperature and high heat environment.
- Do not immerse this product in water or leave it in a damp environment for a long time. Please use a soft cloth to clean the instrument.
- Before use, please confirm that the appearance of this product is not damaged and that there is no abnormality in the boot and screen display. Do not use if there is a malfunction.

Warranty Card

In order to provide a better service, please perfect the following user information:

Within warranty period	Date of purchase:	
Customer Info	Customer Name:	
	Customer address:	
	Tel. No.:	
	Purchase ID:	
	Vendor:	
Warranty Reason		

Warranty Card

1. Please keep this warranty card properly for maintenance certificate.
2. This product is guaranteed for 1 year from the date of purchase.
3. During the warranty period, under normal use and maintenance, the Warranty equipment cannot work because of materials and processes problems. After verification, our company will provide free maintenance and spare parts.
4. During the warranty period, if the following matters occur, the company has the right to refuse service or charge material service fees as appropriate.
 - (1) This warranty card or valid purchase certificate cannot be provided.
 - (2) Product failure and damage caused by improper use by users.
 - (3) Artificial damage caused by abnormal external force.
 - (4) Failure and damage caused by maintenance services not authorized by our company.
 - (5) Failure and damage caused by natural disasters or other force majeure factors.
 - (6) Others intentionally damaged it.
5. The company reserves the right to modify and interpret all content.


Notes



CE certification



FCC certification

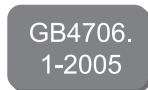

Safety Standard for
Household Appliances



SMART SKIN TECHNOLOGY.

genio
SMART XS

www.genio-innovations.com

© genio Innovations 10/2022 rev. A